

Medical gas emergency kit



Contains:
Medical Devices
IIb Class
CE 0425

Datasheet

The medical gas emergency kit arises from the need to provide hospitals with a ready-to-use kit to cover health emergencies.

In fact, it may happen that you have to power a department in the event of a failure of one of the mains reducers present, or even have to move a patient having the need to connect a small cylinder for oxygen therapy or aerosol therapy.



The kit, which can be supplied either with a comfortable and shockproof portable case (waterproof, dust-proof, resistant to chemicals in compliance with IP67), or with carpentry to be installed in critical areas (painted with Kill Sept anti-battery paint and Easyentry open), have on its inside compact cylinder reducers with hand connection, with a single pressure gauge for an easy and immediate reading of the gas present in the cylinder. Supplied with a low pressure medical gas hose, a medical device clamped with coupling according to Afnor or Uni standards, with a specific gas length of 1.5 meters, which can easily be used to connect the reducers to the gas system, lung ventilators, or any other device that requires FU air and FU Oxygen for its operation.

Portable cases

V-KEM-O2AC-U

STORM emergency kit medical gases for Oxygen and Air FU, NR in a portable shock-proof protector case, including one bottle reducer with Afnor outlet and pressure gauge for Oxygen, one bottle reducer with Afnor outlet and pressure gauge for Air, one 1.5 metre pre-clamped hose with Afnor-UNI coupling for Oxygen (IA-IU), one 1.5 metre pre-clamped hose with Afnor-UNI coupling for Air (IA-IU).RIDT25B+RIDT40B+TUB01015IAIU+TUB0301-SIAIU

V-KEM-02-U

STORM emergency kit medical gases for Oxygen, in portable shock-proof protector case, including one cylinder reducer with one pressure gauge and Afnor outlet for oxygen, one 1.5 metre hose pre-clamped with Afnor-UNI coupling for Oxygen (IA-IU).RIDT25B+TUBE1015IAIU

V-KEM-02AC-A

STORM emergency kit medical gases for Oxygen and Air FU, NR in a portable shock-proof protector case, including one bottle reducer with Afnor outlet and pressure gauge for Oxygen, one bottle reducer with Afnor outlet and pressure gauge for Air, one 1.5 metres pre-clamped with Afnor-Afnor connection for Oxygen (LA-IA), one 1.5 metres pre-clamped hose with Afnor-Afnor connection for Air (IA-IA).RIDT25B+RIDT40B+TUB0301SIAIATA+TUB0101SIAIA

V-KEM-02-A


STORM emergency kit medical gases for Oxygen, in portable shock-proof protector case, including one cylinder reducer with one pressure gauge and Afor oxygen outlet, one 1.5 metre pre-clamped hose with Afnor-Afnor coupling for Oxygen (IA-IA).RIDT25B+TUB01015IA

Outdoor box versions are also available

Technical Data

Typical Dimensions & Weight

- Outside: 41,1 × 32,3 × 16,8 cm
- 5 kg

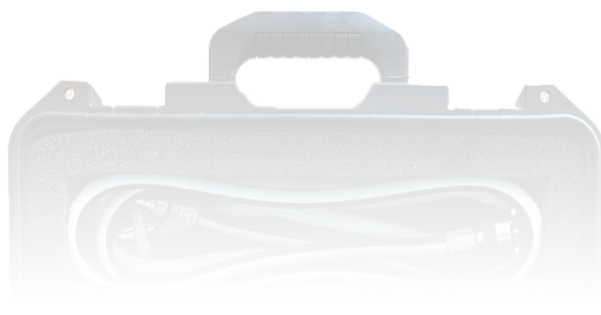



Contacts

DZ Medicale s.r.l

(+39)0307268254

www.dzmedicalesrl.com

info@dzmedicalesrl.com

Via Lega Lombarda 8, 25036,
Palazzolo sull'Oglio (BS) Italia