


Bed Head Unit Girasole

MD


Medical Devices

CE 0425


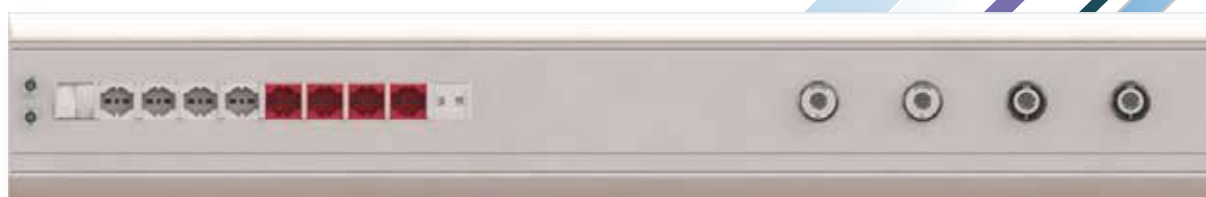


Data sheet

"Girasole" bed head unit, for lighting and medical gas equipment, is made with extruded aluminum profiles and subsequently painted with non-toxic powders.

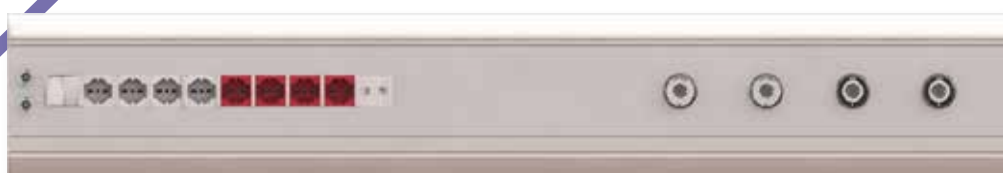



The gas channel, with copper pipes CE Medicalfein EN 13348 strictly isolated from the electrical one, can also be equipped with gas terminal units that follow the standards commonly used today, such as the French one NF S 90-116, to DIN 13260, the UNI 9507 or SS 875 24 30.

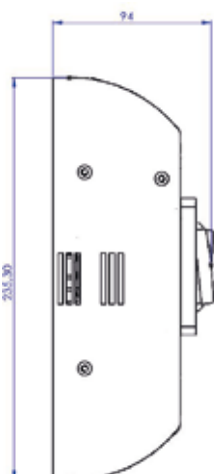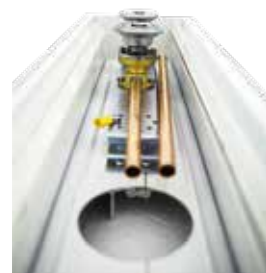

The lighting and electrical composition can be widely customized, and includes ambient lighting, for reading and visiting, and night lighting, a light point suitable for night rest.


built in Norms: EN 7396-1, EN11197, EN 60601-1, EN 60601-1-2

LIGHTING OPTIONS

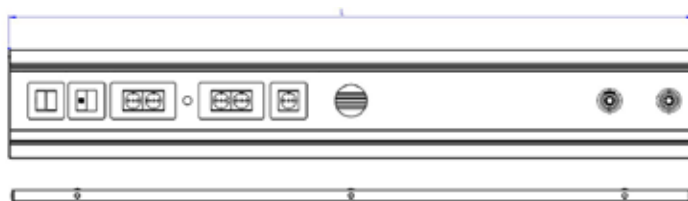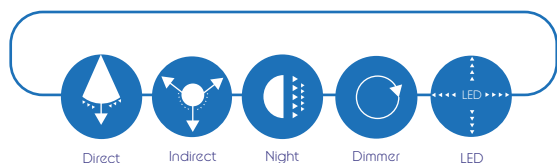

All the LEDs used passed the highest electrical safety standards and the diffuser covering the illuminating part has been chosen with a class of fire resistance greater than required by current standards. The electrical supply points present, which can be built on separate and individually protected lines to keep safe patients and hospital staff, follow different standards such as CEE 7/4 (Schuko type F), UNEL 2P + T 10 / 16A bypass, the French two-pole type E, or BS 1363.


CONFIGURATIONS

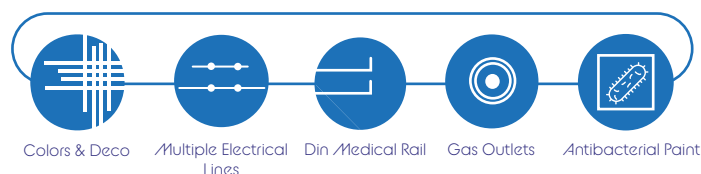


All the presented model can be equipped with a rotating rod, with a hook for phlebo bags and with a standard 30×10 rail that is assembled on the bottom part of the profile. Hand-set (one, two and three pushbutton) for nurse call and lighting control, equipotential socket, pull switch RJ11 and RJ45 data connectors, and other equipment that can be supplied on request. Plates and switches are available also with antibacterial treatment.

Custom Colors


RAL 3012

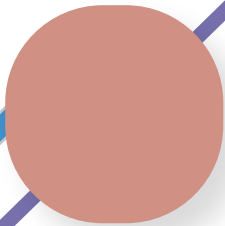


RAL 3011

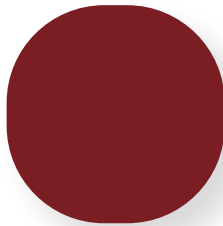


RAL 5020

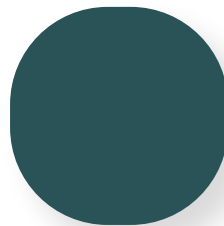


RAL 5012

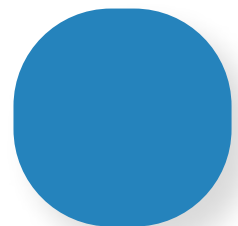


RAL 4008

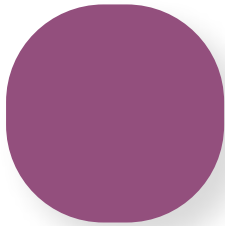


RAL 1016

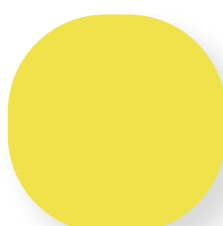


RAL 5024

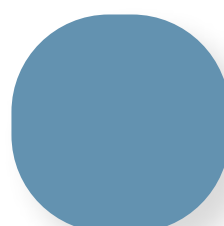


RAL 9006

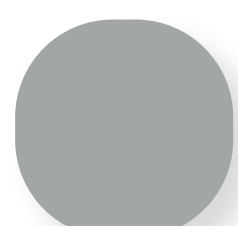


RAL 6019

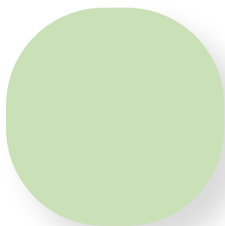


RAL 2000

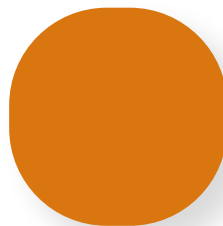


RAL 3020

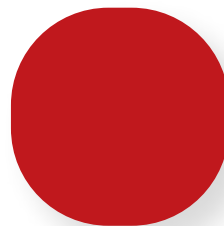


RAL 4009

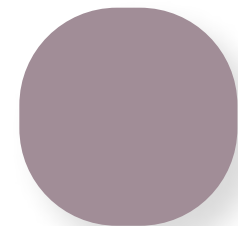


RAL 6018

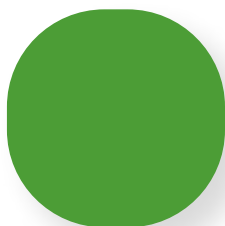

RAL 4003

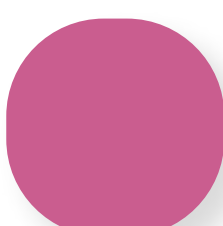

RAL 4003

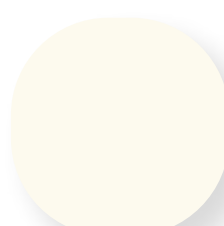

Custom Deco

Adult

City Skyline



Art



Abstract



Maternity



Pediatric

Kids Movies



Comics





Indirizzi utili

DZ Medicale s.r.l

(+39)0307268254

www.dzmedicalesrl.com

info@dzmedicalesrl.com

Via A. De Gasperi n 77, Zocco
di Erbusco (BS) Italia

Via Lega Lombarda 8, 25036,
Palazzolo sull'Oglio (BS) Italia